



User Manual OS 3D PRO DREAMER

Thank you for purchasing the OS 3D PRO DREAMER Massage Chair. Before using your massage chair, please read the contents of this instruction manual and refer to it when needed.

Table of Contents

Features	2
Safety Precautions	3-4
Chair Components	5-6
Chair installing	7-13
Preparation for Use	14
Operation Guide	15-26
Grounding Instructions	27
Chair Specification	28
Troubleshooting & Maintenance	29
Warnings	30

This chapter summarizes the most common problems you could encounter with the massage chair. If you are unable to solve the problem with the information below, contact

If the problem encountered is not listed here, please first try the following steps:

Possible Cause

The appliance is not plugged in.

The appliance has just finished its

automatically switched to standby

15 minutes massage and has

You have used the appliance

continuously for 45 minutes,

to switch off the appliance.

causing the overheat protection

An excessive weight has been put

on the appliance (e.g. someone has

stood on the appliance or a heavy

object has been put on top of it).

You have let the appliance run for

more than 45 minutes in a row.

This sound is caused by the

mechanical structure of the

appliance.

10 Warnings

objects.


Do not use water and other cleaning detergents on the massage chair.
Do not stack, stab, slit or harshly pull on the surface of the massage chair.
Gently press the buttons during operation and avoid hitting the buttons with hard

Pull out the power supply plug after each use.Do not use the massage chair for more than 30 minutes at a time.

• DO NOT position any body part in between the gaps of the leg rest when the automatic leg extension feature is in operation as demonstrated PIC 1, in order to prevent pinch to you. If you are remiss to position any body part in between the gaps of leg rest, please press ON/OFF button to stop and then get your body part out rapidly as demonstrated PIC 2.

• If you have to rise from the chair to stop accident, injury or damage to others by you or you want to get your feet and calves out when the automatic leg extension feature is in operation, please press ON/OFF button to stop as demonstrated PIC 2, the footrest will gradually deflate in some seconds and then quickly get your foot and calves out and rise.

Note: This way as above mentioned just using in emergencies. Please bear in mind that always return the footrest to its original length and restore it to its original position, then rise from the chair in order to prevent accident, injury or damage to yourself and chair.

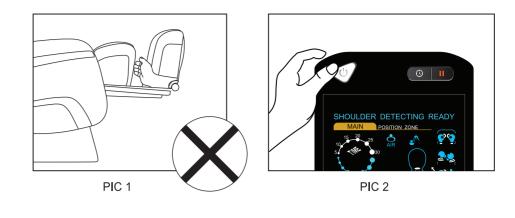

30

2

01 Features

Solution

Put the main plug in the wall socket. Check if the voltage

on the appliance corresponds

If you want another massage,

press the POWER button again.

Unplug the appliance and let it

Unplug the appliance, remove

the excessive weight if it is still

Unplug the appliance and let it

cool down. Always unplug the

This sound is normal and does not mean that the appliance is

appliance after use.

malfunctioning.

on the appliance again and

you use it again.

switch it on.

cool down for 60 minutes before

to the local mains voltage.

The operation of massage chair is controlled by a microcomputer that performs various intelligent massage methods and combinations, like Kneading, Tapping, Clapping,Shiatsu, Rolling, Vibration and Air Massage. All these methods massage acupressure points accurately to effectively alleviate fatigue, promote metabolism and improve blood circulation.

- Equipped with 3D intelligent massage mechanism, user's body curve and acupoints on spine can be positioned more accurately.
- Zero gravity chair.
- Uses 3D rollers to mimic the feel of human thumb.
 34 air bags around the neck, shoulder, back, seat, arms, palms and calves for
- extra comfort as it cradles the body like a glove for increased massage effectiveness.Automatically detect the whole body curve as well as make a micro adjustment
- which brings more scientific massage enjoyment.
- Music sync with two speakers in the backrest.
- New stretching auto program combine with air massage for whole body.Heating therapy in the shoulder.
- Unique feet massage with rollers
- Manual massage for the upper body with three options; Six methods with levels of speed, intensity and three kinds of width adjustable respectively.
 LED light in the side panel create one kind of nice environment.
- The LCD displayer.
- Auto-timer 5-30 minutes options.

29

1

09

immediately.

Troubleshooting

1. Unplug the product for 60 minutes.

3. Turn on the Power Switch.

Problem

The massage chair

does not work.

The appliance

Sound can be heard

during massage.

is warm.

2. After 60 minutes, plug in the product again

mode.

02 Safety Precautions

the gaps of the massage chair.

When using the massage chair, basic precautions should always be followed to reduce the risk of electric shock, burns, fire or injury to persons. Please read the following precautions before using the massage chair.Make sure the power plug is securely inserted into a proper electrical outlet to

reduce the risk of short-circuit and fire.
Always unplug the massage chair from the electrical outlet immediately after using and before cleaning.
DO NOT use the massage chair in areas where the humidity is high such as a

DO NOT use the massage chair outdoors and avoid areas with direct sunlight.

If the electricity is off, turn the power switch to the OFF position and unplug the massage chair from the electrical outlet to prevent any damages to the chair
when the electricity comes back on.

Never plug/unplug the power cord plug socket with wet hands.
Close and continuous supervision is necessary when the massage chair is used by, on or near children or disabled persons.

No more than one person should use the massage chair at a time.
The massage chair is intended for use as described in the Instruction Manual. Do not use any attachments other than those recommended by the manufacturer.
DO NOT attempt to stand or sit on the backrest, armrest or footrest in order to

prevent accident, injury or damage to the chair.
Avoid falling asleep while using the massage chair.
Always keep the massage chair clean, and never insert any objects into or around


DO NOT move the massage chair by pulling on the power supply cord.
Make sure there are no children, pets or other obstacles behind or underneath the massage chair.

• If you find any rips or tears in the fabric of the massage chair contact the service department and stop using the chair immediately.

08 Product Specifications

Description	Specifications		
Model No.	OS 3D PRO DREAMER		
Rated Voltage	110-120V~		
Rated Frequency	50/60Hz		
Rated Power Consumption	250W		
Timer	5/10/15/25/30 minutes		
Import sound source	MP3, MP4, mobile phone		
The power of loudspeaker	2W (one side)		
Dimension (L x W x H)	Backrest : 54.1" x 30.1" x 23.0"		
	Chair Base: 50.2" x 30.3" x 25.6"		
	Side panel : 25.8" x 15.7" x 39.4"		
Weight	Product gross weight: 306. 4lbs		
Weight	Product net weight: 261.2 lbs		
	Controller wire: 59"		
Length of the wire	Power supply wire: 71"		
	Audio wire: 25.6"		
Llagra condition	Environment temperture: 10°-40°		
Usage condition	Contrasting humidity: 30-85RH		
Otana and dition	Storage temperture: 20°-60°		
Storage condition	Storage humidity level: 30-85RH		
Safety feature	Equipped with overheated and power surge safety protection.		
Using benefits	Increasing blood circulation; Relieving muscular fatigue		



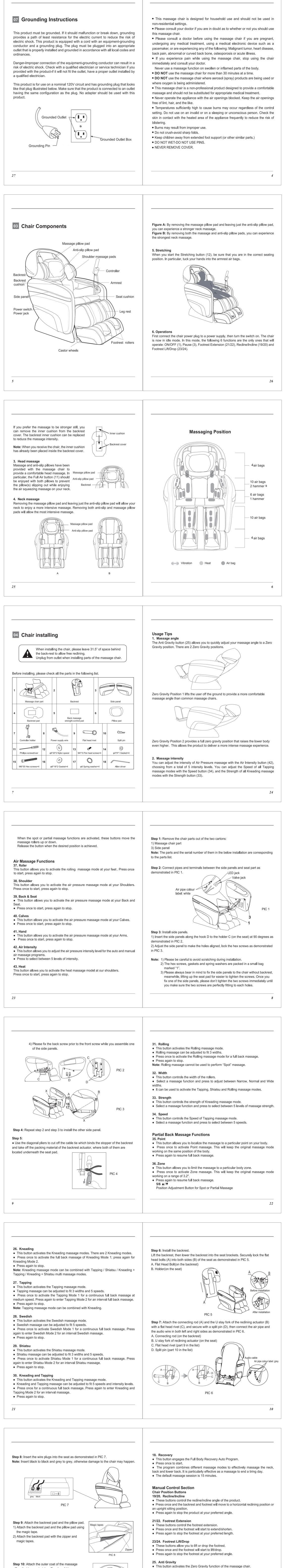

Step 10: Attach the outer coat of the massage pillow that is easy to take off and could be cleaned by water.

PIC 9

11

20

• A Zero Gravity position is a recline position which mimics the posture of astronauts


during lift off on space missions. This position is designed to distribute weight effectively across the chair allowing the maximum in stress relief. This also allows for a more intense massage experience. For some users, a Zero Gravity position allows them to feel

• Press once to activate Zero Gravity Position 1. Press again to activate Zero Gravity

weightlessness.

Press again to stop.

Position 2.

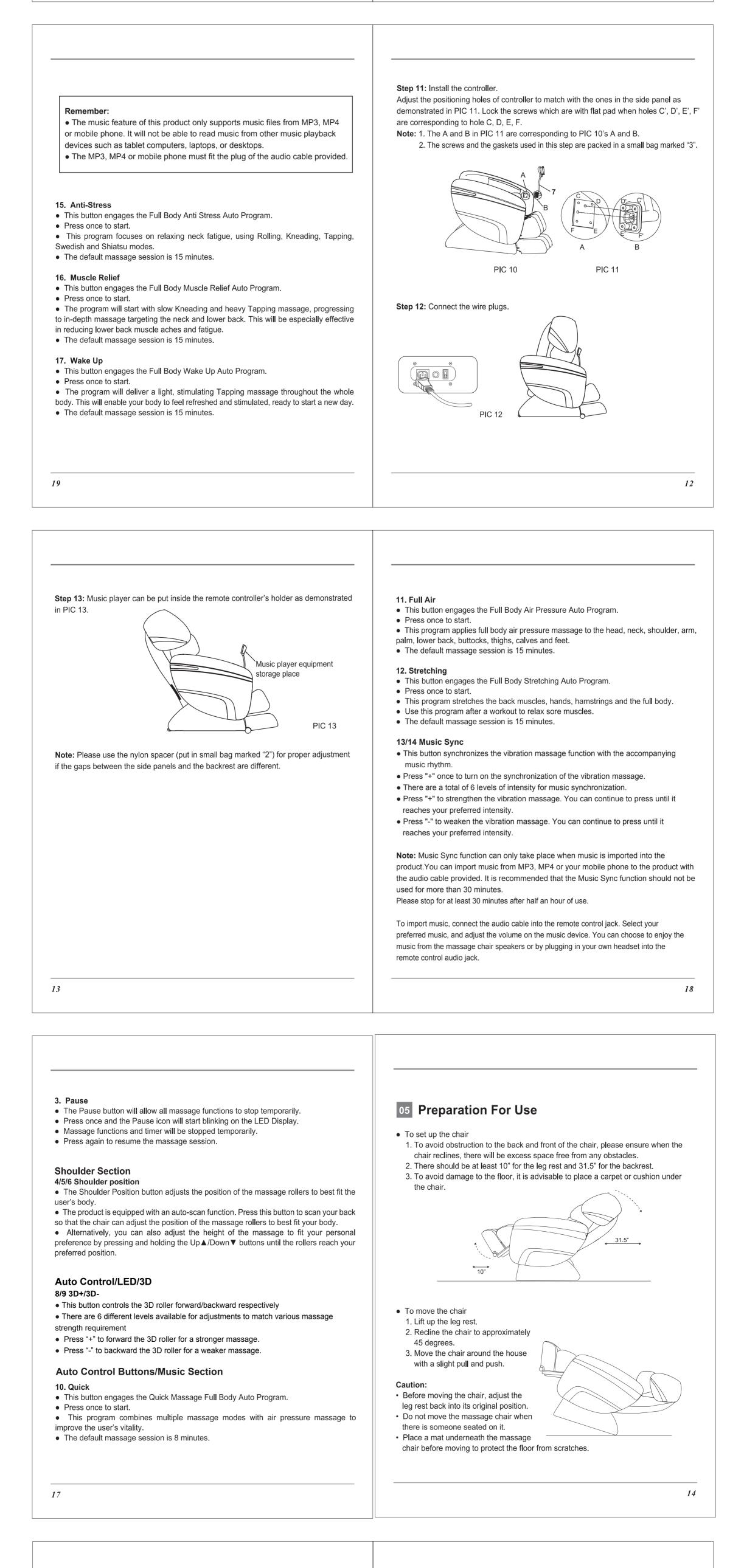



The Remote Control is divided into 4 main sections with different functional purposes. The pages that follow describe the functions of the various buttons within these sections. **Basic Control Section**

Section	Buttons	Reference	Section	Buttons	Reference
Basic Control	Power	1		Chair Position	19-24
	Timer	2		Zero Gravity Position	25
		_		Massage Mode	26-31
	Pause	3		Massage Strength	33-34
Shoulder	Shoulder Position	4-6	Manual	Width	32
Auto Control / LED / Music	LED	7	Control	Point / Zone Massage	35-36
	Music Sync	8-9		Air Massage	37-41
	Auto Programs	10-18		Air Intensity	42

Connect to the power supply, and turn the switch to "I" position. The chair is now in idle mode. During idle mode, there are only seven functions: ON/OFF, Foot Extension, Pause, Recline, Incline, Footrest lift, Footrest drop that can be operated.

Basic Control Section 1. Power

- The Power button starts and stops all functions on the massage chair.Press once.
- Remote control and side panel LED lights will start blinking to indicate the chair is in idle mode.
- You can press any massage function buttons now to begin your massage session, at a default time of 15 minutes.
- Press again. All functions will stop, and massage rollers will return to original position.
- 2. Timer
- The Timer button controls the length of the massage session.
 The default time for all massage functions is 15 minutes.
- Press each time to change the length of time in the following sequence: 20>25>30>05>10>15 minutes.